
HAMBLINS UNIVERSAL SLIT-LAMP (No. 620) AND TABLE (No. 602/c)

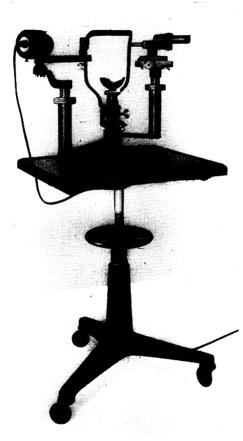

No. 620

Slit Lamp complete with Corneal Microscope with one pair of objectives, head and chin rests.

No. 602. Table (without castors). No. 602/c. Table (with castors).


In Hamblin's Universal Slit-Lamp and Corneal Microscope the principle of having both the microscope and the lamp each borne on independent arms swinging from common central upright support was employed from the outset. This principle now generally adopted in instruments of this kind makes for portability and allows the apparatus to be used on any table adjustable for height.

In the last quarter of a century the Universal Slit-Lamp has become notable for ease in manipulation, accuracy and smoothness of movement, and for intrinsic workmanship and lasting quality.

It can be fitted if required at a small extra charge with the Goldmann Angle Reduction Prism for easier examination of the posterior parts of the eye.

A suitable table is shewn in the illustration which is quickly adjusted for height, easily propelled and solidly built.

THEODORE HAMBLIN L^{TD} DISPENSING OPTICIANS 15 WIGMORE STREET, I ONDON,W.1.

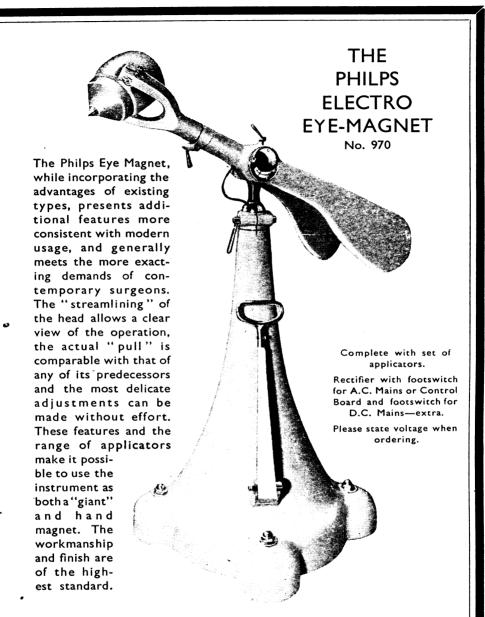

Operating Theatre regulations require flameproof switch gear to be used with all apparatus of this type. Where A.C. only is available a Rectifier is necessary. This can be supplied complete with suitable footswitch at extra cost.

THEODORE
HAMBLIN LTD
DISPENSING OPTICIANS
15 WIGMORE STREET
IONDON.W.1

GEORGE SPILLER LTD.

dispensing opticians
WIGMORE STREET, W.I.

'Phone: Welbeck 7471-2.

EVEN TINT.

Vacuum deposited "Bloom"

The extension of this well known process to sighted Ophthalmic lenses, provides the ideal solution of the problem of protection against excess light.

It is available in several degrees of absorption and is even over the lens whatever the prescription.

